




# PROJECTION **MICROSCOPE**

MODEL NO. - ACM- 22306-J



# PROJECTION **MICROSCOPE**

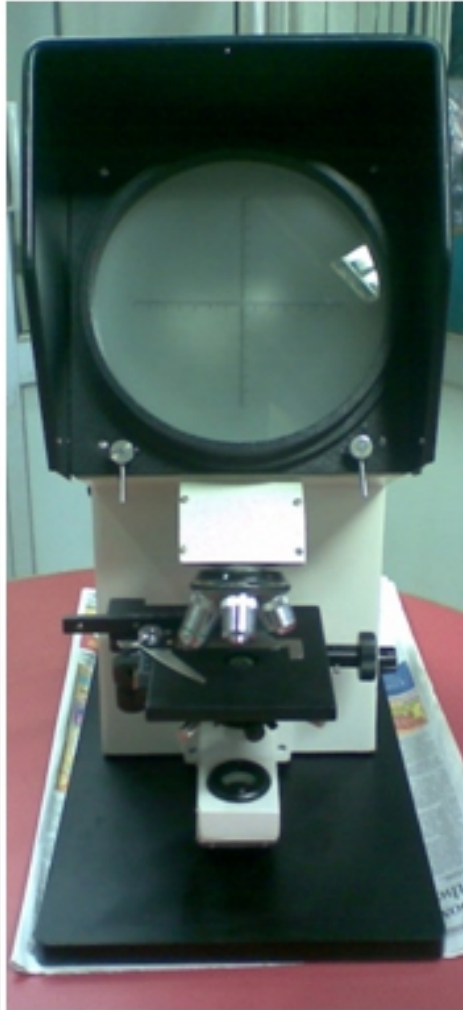

## **SPECIFICATIONS**

An advanced versatile projection microscope for measurement of specimens, different types of yarns in textile industries and research institutes. It is equipped with Halogen bulb of 12V/100 watt workable through a solid state built in transformer working on 220V AC mains. 200mm. dia graduated screen rotatable through 360°, built in fresnel lens with loop and arrow arrangement, coaxial mechanical stage. Built in adjustable compensating eyepiece is connected with magnification correction knob. Focussing through coarse and fine motion knobs. Quadruple nose-piece with four objectives-5x, 10x, 20x, 40x. Supplied in sturdy cardboard box wt. approx. 20kg.



An ISO 9001 : 2208 | ISO 14001 : 2008 | ISO 13485  
WHO: GMP Products | GLP Compliant Products

[www.acmasindia.com](http://www.acmasindia.com) | [www.acmasglobal.com](http://www.acmasglobal.com) | [www.test-chambers.com](http://www.test-chambers.com)

## **ACMAS TECHNOLOGIES PVT. LTD.**

### **CORPORATE OFFICE/ SHOWROOM (INDIA)**

Plot No. 352-353, Sector – 57  
Phase- IV Kundli, Sonapat, Haryana - 131028  
Land Line No.: +91 - 11- 47619688,  
Fax: +91-11-47619788  
E-mail.: [info@acmasindia.com](mailto:info@acmasindia.com)

### **SALES OFFICE (HONG KONG)**

Unit D 28 11/F Wing Tat Comm,Bidg 97, Bonham Strand East,  
Sheung Wan, Hong Kong (PRC)  
Tel.: 0086-13929598046 | 0086- 18922303099  
E-mail.: [hk@acmasindia.com](mailto:hk@acmasindia.com)

### **SALES OFFICE (RUSSIA)**

Inmed Trade Street Ozerkovsky Embankment,  
Unit No 50, Straine- 1, Off- 502, Moscow, Russia  
E-mail.: [info@acmasindia.com](mailto:info@acmasindia.com)  
Tel.: 0049- 79592345 | Email: [russia@acmasindia.com](mailto:russia@acmasindia.com)